


SEQUENCE FOR **PUTTING ON** PERSONAL PROTECTIVE EQUIPMENT (PPE) with CAPR for COVID AREA

Make sure to always use proper hand hygiene before each step


Training

- Log into testbank and click “select test”
- Under training type: Clinical Training, select “CAPR Training for non-OR staff (1692)”
- Obtain CAPR Lens Cover (DLC) and Comfort band from manager/director once verified

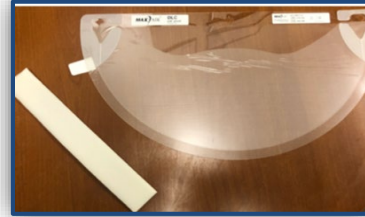

Battery Pack/Belt

- Obtain a fully charged battery (LED should be green on charger)
- Assemble battery to belt using clip
- Place belt comfortably around waist, with battery near side/back of hip

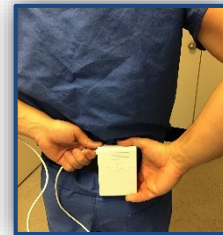


Lens and Battery Cord

- Attach foam comfort strip and DLC lens to helmet, with DLC cuff facing the inside of the helmet
- Inside helmet, pull plastic temple flapper away from lens to rest on sponge material
- Remove inner and outer film if not already removed
- Connect power cord to battery pack, listen for audible click of connection
- 5 lights will illuminate,
 - Green indicator should remain on
 - Amber indicator unit needs inspected
 - Red indicator = battery 25% or less



Foam comfort strip

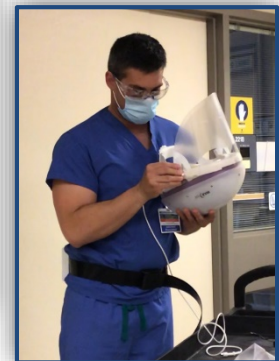

Donning CAPR Unit and PPE

- Loosen adjustment knob so helmet easily fits over head
- Hold helmet in one hand, while pulling down the cuff under chin, and place helmet on head
- Assure good seal by sliding fingers up and down between face and cuff.
 - Adjust knob until helmet is stable on head
- Put on hospital approved gown and gloves
 - Gown should cover CAPR battery/cord



Always use proper hand hygiene before each step





SEQUENCE FOR REMOVING PERSONAL PROTECTIVE EQUIPMENT (PPE) with CAPR for COVID AREA

ALWAYS USE PROPER HAND HYGIENE WHILE REMOVING PPE


Gown and Glove Removal

- Remove gloves and gown inside of patient room
- Perform hand hygiene
(Alcohol based hand rub or soap and water)

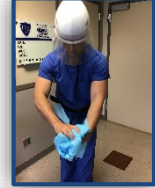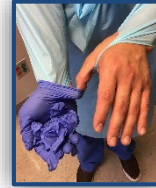


CAPR Lens

- Proceed into hallway to CAPR cart (near MICU door) to doff CAPR
- Perform hand hygiene
- Unplug cord from battery pack
- Perform hand hygiene
- Remove DLC lens from CAPR system
- Use disinfectant wipe to clean the entirety of the DLC lens.
- Perform hand hygiene
- Place lens into manila envelope with Name and D.O.B.

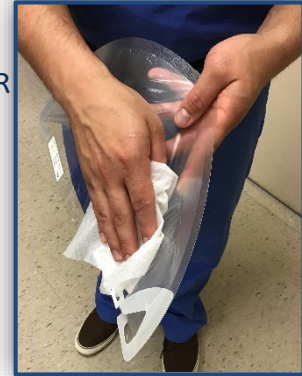

CAPR Battery, Belt, Helmet

- Perform hand hygiene
- Clean the battery pack and belt
 - Use hospital approved disinfectant wipe
 - Place onto CAPR supply cart
- Perform hand hygiene
- Remove comfort band strip and set on manila folder
- Remove CAPR and use hospital approved disinfectant wipe
 - Start from the inside and work towards the outer shell
- Grab new wipe if needed for proper cleaning/disinfecting

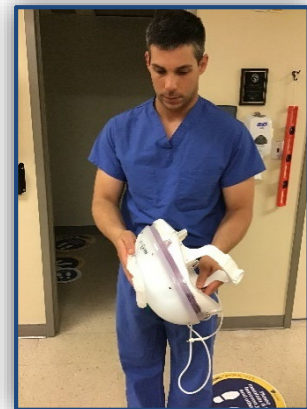

CAPR Helmet/Battery Pack

- Perform hand hygiene
- Return CAPR system back to cart
- Reconnect battery pack to charging station
- Perform hand hygiene
- Exit space and put on proper PPE (mask/eye protection)



Always use proper hand hygiene before each step

